

Bausch & Lomb SofLens®

daily disposable

Toric (hilafilcon B)

For Astigmatism

Visibility Tinted Contact Lenses

CAUTION: Federal (U.S.A.) law restricts this device to sale by or on the order of a licensed practitioner.

Bausch & Lomb Incorporated
Rochester, NY 14609

© Bausch & Lomb Incorporated. All rights reserved worldwide.
®/TM are trademarks of Bausch & Lomb Incorporated or its affiliates. Printed in U.S.A.

SL7496
8077101

DESCRIPTION

The Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens is a soft hydrophilic contact lens which is available in a flexible shell with a toric surface. The lens is made from the hilafilcon B material, a hydrophilic copolymer of 2-hydroxyethyl methacrylate and N-vinyl pyrrolidone, and is 59% water by weight when immersed in a sterile saline solution. This lens is tinted blue with Reactive Blue Dye 246.

The physical / optical properties of the lens are:

Specific Gravity:	1.119
Refractive Index:	1.4036
Light Transmittance:	C.I.E. Y value - approximately 95%
Water Content:	59%
Oxygen Permeability:	$22 \times 10^{-11} \text{cm}^3 \text{O}_2 (\text{STP}) \times \text{cm} / (\text{sec} \times \text{cm}^2 \times \text{mmHg}) @ 35^\circ \text{C}$ (Polarographic Method)


The Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens is to be prescribed for single-use disposable wear.

The Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lenses are manufactured by a cast molding process. The anterior surface of the lens contains the spherical power, prism ballast and spherocylindrical curve to accommodate the required astigmatic power.

SYMBOL REFERENCE GUIDE

For label and cartons:

	Quality System Certification Symbol		Caution: Federal (U.S.A.) law restricts this device to sale by or on the order of a licensed practitioner
	Fee Paid for Waste Management		Single Use Only
	Sterile Using Steam or Dry Heat		Base Curve
	See Instruction Leaflet		Sphere Power (Diopters)
	Diameter		Cylinder Axis (Degrees)
	Use by Date (Expiration Date)		Cylinder Power (Diopters)
	Diopter (Lens Power)		Batch Code
	Authorized Representative in European Community		

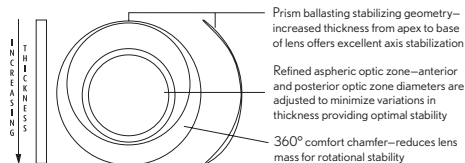

Guide Mark System

Each Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens is marked with a Guide Mark in the lens perimeter at 6 o'clock. This Guide Mark gives an instant reference for estimating lens rotation and orientation. It is, in effect, a protractor guide on the lens surface. The guide marking makes proper axis orientation and fitting faster and easier. The guide marking need not position at exactly 6 o'clock on the cornea.

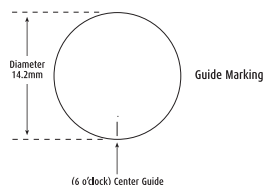


TABLE OF CONTENTS

Introduction	2
Description	3
Lens Parameters Available	5
How the Lens Works (Actions)	5
Indications	5
Contraindications (Reasons Not To Use)	6
Warnings	7
Precautions	8
Adverse Reactions	10
Important Treatment Information for Adverse Reactions	11
Selection of Patients	12
Fitting Procedure	13
Pre-Fitting Examination	13
Initial Lens Power Selection	13
Initial Lens Evaluation	14
Criteria of a Well-Fitted Lens	15
Characteristics of a Tight (Steep) Lens	15
Characteristics of a Loose (Flat) Lens	15
Follow-up Care	16
Professional Fitting Sets	17
Wearing Schedule	17
Monovision Fitting Guidelines	17
Patient Selection	17
Eye Selection	18
Special Fitting Considerations	19
Near Add Determination	19
Trial Lens Fitting	19
Adaptation	20
Other Suggestions	20
Handling of Lenses	21
Care for a Sticking (Nonmoving) Lens	21
Emergencies	21
Reporting of Adverse Reactions	22
How Supplied	22

LENS PARAMETERS AVAILABLE

The Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens is a hemispherical shell of the following dimensions:

Diameter:	14.2mm
Center Thickness:	0.05mm to 0.75mm
Base Curve:	8.6mm
Powers (Spherical):	plano to -6.00D in 0.25D steps -6.50D to -9.00D in 0.50D steps
Cylinder Powers:	-0.75D, -1.25D, and -1.75D
Axis:	0° to 180° in 10° Increments

HOW THE LENS WORKS (ACTIONS)

In its hydrated state, the Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens has a unique Lo-Torque™ design that results in excellent stability, visual acuity and comfort. When placed on the cornea, the Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens acts as a refracting medium to focus light rays on the retina.

INDICATIONS

The Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens is indicated for the daily wear correction of refractive ametropia (myopia, hyperopia, and astigmatism) in not-aphakic persons with non-diseased eyes, exhibiting astigmatism of 5.00 diopters or less, that does not interfere with visual acuity. The lens may be prescribed in spherical powers ranging from +20.00D to -20.00D.

The lens is to be prescribed for single-use disposable wear, and is to be discarded after each removal.

INTRODUCTION

This package insert has been developed to provide eye care professionals with information covering characteristics of the Bausch + Lomb SofLens® daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens and to illustrate fitting procedures. It is effective as of February 2015 and supersedes all prior fitting guides for the product described. Please read carefully and keep this information for future use.

This package insert is intended for the eye care professional, but should be made available to patients upon request. The eye care professional should provide the patient with the patient instructions that pertain to the patient's prescribed lens, and the recommended wearing schedule.

CONTRAINDICATIONS (REASONS NOT TO USE)

DO NOT USE the Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens when any of the following conditions exist:

- Acute and subacute inflammation or infection of the anterior chamber of the eye
- Any eye disease, injury, or abnormality that affects the cornea, conjunctiva, or eyelids
- Severe insufficiency of lacrimal secretion (dry eyes)
- Corneal hypoesthesia (reduced corneal sensitivity)
- Any systemic disease that may affect the eye or be exaggerated by wearing contact lenses
- Allergic reactions of ocular surfaces or adnexa (surrounding tissue) that may be induced or exaggerated by wearing contact lenses or use of contact lens solutions
- Any active corneal infection (bacterial, fungal, or viral)
- If eyes become red or irritated

WARNINGS

- After a thorough eye examination, including appropriate medical background, patients should be fully apprised by the prescribing professional of all the risks with contact lens wear. Patients should be advised of the following warnings pertaining to contact lens wear:
 - Problems with contact lenses could result in **serious injury** to the eye. It is essential that patients follow their eye care professional's direction and all labeling instructions for proper use of lenses. Eye problems, including corneal ulcers, can develop rapidly and lead to **loss of vision**.
 - Daily wear lenses are not indicated for overnight wear, **and patients should be instructed not to wear lenses while sleeping**. Clinical studies have shown that the risk of serious adverse reactions is increased when daily wear lenses are worn overnight.
 - Studies have shown that contact lens wearers who are smokers have a higher incidence of adverse reactions than nonsmokers
 - If a patient experiences eye discomfort, excessive tearing, vision changes, or redness of the eye, the patient should be instructed to **immediately remove lenses** and promptly contact his or her eye care professional.

7

IMPORTANT TREATMENT INFORMATION FOR ADVERSE REACTIONS

Sight-threatening ocular complications associated with contact lens wear can develop rapidly, and therefore early recognition and treatment of problems are critical. Infectious corneal ulceration is one of the most serious potential complications, and may be ambiguous in its early stage. Signs and symptoms of infectious corneal ulceration include discomfort, pain, inflammation, purulent discharge, sensitivity to light, cells and flare and corneal infiltrates.

Initial symptoms of a minor abrasion and an early infected ulcer are sometimes similar. Accordingly, such epithelial defect, if not treated properly, may develop into an infected ulcer. In order to prevent serious progression of these conditions, a patient presenting symptoms of abrasions or early ulcers should be evaluated as a potential medical emergency, treated accordingly, and be referred to a corneal specialist when appropriate. Standard therapy for corneal abrasions such as eye patching or the use of steroids or steroid/antibiotic combinations may exacerbate the condition. If the patient is wearing a contact lens on the affected eye when examined, the lens should be removed immediately and the lens and lens care products retained for analysis and culturing.

11

PRECAUTIONS

Precautions for Eye Care Professionals

- Due to the small number of patients enrolled in clinical investigation of lenses, all refractive powers, design configurations, or lens parameters available in the lens material are not evaluated in significant numbers. Consequently, when selecting an appropriate lens design and parameters, the eye care professional should consider all characteristics of the lens that can affect lens performance and ocular health, including oxygen permeability, wettability, central and peripheral thickness, and optic zone diameter.

The potential impact of these factors on the patient's ocular health should be carefully weighed against the patient's need for refractive correction; therefore, the continuing ocular health of the patient and lens performance on eye should be carefully monitored by the prescribing eye care professional.
- Eye care professionals should instruct the patient to **REMOVE A LENS IMMEDIATELY** if an eye becomes red or irritated.
- Fluorescein, a yellow dye, should not be used while the lenses are on the eyes. The lenses absorb this dye and become discolored. Whenever fluorescein is used in the eyes, the eyes should be flushed with sterile saline solution that is recommended for in-eye use.
- The patient should be instructed to always discard disposable lenses.
- Aphakic patients should not be fitted with Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lenses until the determination is made that the eye has healed completely.
- The lenses are prescribed for disposable wear, and are to be disposed of once they are removed from the patient's eye. It is important that patients be instructed to always have available a pair of replacement lenses. In the event that a lens must be removed from the eye because of dust, a foreign body or other contaminant gets on the lens or the lens becomes dehydrated, the lens should be removed and replaced with a replacement lens.
- Eye care professionals should carefully instruct patients about the following safety precautions. It is strongly recommended that patients be provided with a copy of the Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens Patient Information Booklet available from Bausch + Lomb and understand its contents prior to dispensing the lenses.

8

SELECTION OF PATIENTS

The eye care professional should not fit patients who cannot or will not adhere to the recommended replacement regimen or are unable to place and remove the lenses should not be provided with them. Failure to follow handling instructions could lead to serious eye infections which might result in corneal ulcers.

Patient communication is vital because it relates not only to patient selection but also to ensure compliance. It is also necessary to discuss the information contained in the Patient Information Booklet with the patient at the time of the initial examination.

Patients selected to wear Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lenses should be chosen for their motivation to wear contact lenses, general health and cooperation. The eye care professional must take care in selecting, examining and instructing contact lens patients. Patient hygiene and willingness to follow the eye care professional instructions are essential to their success.

A detailed history is crucial to determining patient needs and expectations. Your patient should be questioned regarding vocation, desired lens wearing time (full or part time), and desired lens usage (reading, recreation or hobbies).

Initial evaluation of the lens should be preceded by a complete eye examination, including visual acuity with and without correction at both distance and near, keratometry and slit lamp examination.

It is normal for the patient to experience mild symptoms such as lens awareness, variable vision, occasional tearing (watery eyes) and slight redness during the adaptation period. Although the adaptation period varies for each individual, generally within one week these symptoms will disappear. If these symptoms persist, the patient should be instructed to contact his or her eye care professional.

12

Handling Precautions

- Always wash and rinse hands before handling lenses. Do not get cosmetics, lotions, soaps, creams, deodorants, or sprays in the eyes or on the lenses. It is best to put on lenses before putting on makeup. Water-base cosmetics are less likely to damage lenses than oil-base products.
- Be sure that before leaving the eye care professional's office, the patient is able to remove lenses promptly or have someone else available to remove them.
- Be certain that the fingers or hands are free of foreign materials before touching lenses, as microscopic scratches of the lenses may occur, causing distorted vision and/or injury to the eye.
- Always handle lenses carefully and avoid dropping them.
- Do not touch the lens with fingernails.
- Carefully follow the handling, insertion, removal, cleaning/disinfecting, storing and wearing instructions in the Patient Information Booklet for the Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lenses and those prescribed by the eye care professional.
- Never use tweezers or other tools to remove lenses from the lens container unless specifically indicated for that use. Pour the lens into the hand.

Topics to Discuss with the Patient

Who Should Know That the Patient is Wearing Contact Lenses:

As with any contact lens, follow-up visits are necessary to assure the continuing health of the patient's eyes. The patient should be instructed as to a recommended follow-up schedule.

- Patients should inform their doctor (health care professional) about being a contact lens wearer.
- Patients should always inform their employer of being a contact lens wearer. Some jobs may require the use of eye protection equipment or may require that you do not wear lenses.

9

FITTING PROCEDURE

1. Pre-Fitting Examination

A pre-fitting patient history and examination are necessary to:

- determine whether a patient is a suitable candidate for daily wear contact lenses (consider patient hygiene and mental and physical state),
- make ocular measurements for initial contact lens parameter selection, and
- collect and record baseline clinical information to which post-fitting examination results can be compared.

A prefitting examination should include spherocylinder refraction and VA, keratometry, and biomicroscopic examination.

2. Initial Lens Power Selection

- a. Select the initial trial lens from the Toric Diagnostic Lens Set with a power most similar to the patient's refractive needs, or order a diagnostic lens to the prescription which most closely matches that of the patient.
- b. Place the lens on the eye and allow the lens to remain on the eye long enough (10 to 20 minutes) to achieve a state of equilibrium. Small variations in the tonicity, pH of the lens solutions, and individual tear composition may cause slight changes in fitting characteristics.
- c. Allow any increase in tear flow to subside before evaluating the lens. The time required will vary with the individual.

13

ADVERSE REACTIONS

The patient should be informed that the following problems may occur:

- Eyes stinging, burning, itching (irritation), or other eye pain
- Comfort is less than when lens was first placed on eye
- Abnormal feeling of something in the eye (foreign body, scratched area)
- Excessive watering (tearing) of the eyes
- Unusual eye secretions
- Redness of the eyes
- Reduced sharpness of vision (poor visual acuity)
- Blurred vision, rainbows, or halos around objects
- Sensitivity to light (photophobia)
- Dry eyes

If the patient notices any of the above, he or she should be instructed to:

- **Immediately remove the lenses.**
- If the discomfort or problem stops, then look closely at the lens. If the lens is in any way damaged, **do not** put the lens back on the eye. Place the lens in the storage case and contact their eye care professional. If the lens has dirt, an eyelash, or other foreign body on it, or the problem stops and the lens appears undamaged, the patient should insert a new lens. After insertion of a new lens, if the problem continues, the patient should **immediately remove the lenses and consult their eye care professional.**

If the above symptoms continue after removal of the lens, or upon insertion of a new lens, the patient should **immediately remove the lenses and contact his or her eye care professional** or physician, who must determine the need for examination, treatment or referral without delay. (See Important Treatment Information for Adverse Reactions.) A serious condition such as infection, corneal ulcer, corneal vascularization, or iritis may be present, and may progress rapidly. Less serious reactions such as abrasions, epithelial stinging or bacterial conjunctivitis must be managed and treated carefully to avoid more serious complications.

10

3. Initial Lens Evaluation

- a. To determine proper lens parameters, observe the lens relationship to the eye using a slit lamp. The toric diagnostic lens is used to optimize lens fitting characteristics and determine axis orientation. Lens power is determined by the spectacle refraction.
 - Rotation evaluation: The center guide mark should locate at the inferior limbus. Once oriented, rotational rocking should be limited to less than 5°.
 - Movement: The lens should provide discernible movement with:
 - Primary gaze blink
 - Upgaze blink
 - Upgaze lag
 - Centration: The lens should provide full corneal coverage
- b. Determine contact lens power. When the toric diagnostic lens does not have a power equivalent to their spectacle Rx, Sphero-cylinder over-refractions will often be inaccurate and confusing. Therefore it is usually preferable to use the spectacle Rx as the only basis for the contact lens power. The sphere and cylinder power of the spectacle Rx becomes the sphere and cylinder power of the contact lens. There are two exceptions:
 1. If spectacle cylinder power falls between available contact lens cylinder powers, prescribe the lesser contact lens cylinder power. The sphere power can be increased -0.25D to compensate if desired. Of course this can vary depending on your interpretation of the patient's subjective responses.

Example:
Spectacle Rx -2.00 -1.00 X 180
Contact Lens Power Ordered -2.25 -0.75 X 180
 2. When the spectacle lens power in any principle meridian is greater than 4.00D, the spectacle refraction should be vertexed to the corneal plane. This can affect both the sphere and cylinder powers ordered.

Example:
Spectacle Rx -5.00 -1.75 X 180
Contact Lens Power Ordered -4.75 -1.25 X 180

14

4. Criteria of a Well-Fitted Lens

If the initial lens selection fully covers the cornea, provides discernible movement after a blink, is comfortable for the patient and provides satisfactory visual performance, it is a well fitted lens and can be dispensed.

5. Characteristics of a Tight (Steep) Lens

A lens which is much too steep may subjectively and objectively cause distortion which will vary after a blink. However, if a lens is only marginally steep, the initial subjective and objective vision and comfort findings may be quite good. A marginally steep lens may be differentiated from a properly fitted lens by having the patient gaze upward. A properly fitted lens will tend to slide downward approximately 0.5mm while a steep lens will remain relatively stable in relationship to the cornea, particularly with the blink.

With your finger, gently rotate the lens approximately 45° to the temporal side. It should reorient within 5 to 10 blinks back to the same stabilized position.

6. Characteristics of a Loose (Flat) Lens

If the lens is too flat, it will:

- Decenter, especially on post-blink.
- Have a tendency to edge lift inferiorly and sit on the lower lid, rather than positioning between the sclera and palpebral conjunctiva.
- Have a tendency to be uncomfortable and irritating with fluctuating vision.
- Have a tendency to drop or lag greater than 2.0mm on upgaze post-blink.

15

3. Special Fitting Considerations

Unilateral Lens Correction

There are circumstances where only one contact lens is required. As an example, an emmetropic patient would only require a near lens while a bilateral myope may require only a distance lens.

Example: A presbyopic emmetropic patient who requires a +1.75 diopter add would have a +1.75 lens on the near eye and the other eye left without a lens.

A presbyopic patient requiring a +1.50 diopter add who is -2.50 diopters myopic in the right eye and -1.50 diopters myopic in the left eye may have the right eye corrected for distance and the left uncorrected for near.

4. Near Add Determination

Always prescribe the lens power for the near eye that provides optimal near acuity at the midpoint of the patient's habitual reading distance. However, when more than one power provides optimal reading performance, prescribe the least plus (most minus) of the powers.

5. Trial Lens Fitting

A trial fitting is performed in the office to allow the patient to experience monovision correction. Lenses are fit according to the directions in the general fitting guidelines.

Case history and standard clinical evaluation procedure should be used to determine the prognosis. Determine which eye is to be corrected for distance and which eye is to be corrected for near. Next determine the near add. With trial lenses of the proper power in place observe the reaction to this mode of correction.

Immediately after the correct power lenses are in place, walk across the room and have the patient look at you. Assess the patient's reaction to distance vision under these circumstances. Then have the patient look at familiar near objects such as a watch face or fingernails. Again assess the reaction. As the patient continues to look around room at both near and distance objects, observe the reactions. Only after these vision tasks are completed should the patient be asked to read print. Evaluate the patient's reaction to large print (e.g. typewritten copy) at first and then graduate to news print and finally smaller type sizes.

After the patient's performance under the above conditions are completed, tests of visual acuity and reading ability under conditions of moderately dim illumination should be attempted.

An initial unfavorable response in the office, while indicative of a guarded prognosis, should not immediately rule out a more extensive trial under the usual conditions in which a patient functions.

19

7. Follow-up Care

a. Follow-up examinations are necessary to ensure continued successful contact lens wear. From the day of dispensing, the following schedule is a suggested guideline for follow up.

- 24 hours
- 10 days
- 1 month
- 3 months
- every six months thereafter

At the initial follow-up evaluations the eye care professional should again reassure the patient that any of the previously described adaptive symptoms are normal, and that the adaptation period should be relatively brief.

b. Prior to a follow-up examination, the contact lenses should be worn for at least 4 continuous hours and the patient should be asked to identify any problems which might be occurring related to contact lens wear.

c. With lenses in place on the eyes, evaluate fitting performance to assure that CRITERIA OF A WELL FITTED LENS continue to be satisfied. Examine the lenses closely for surface deposition and/or damage.

d. After the lens removal, instill sodium fluorescein (unless contraindicated) into the eyes and conduct a thorough biomicroscopy examination.

1. The presence of vertical corneal striae in the posterior central cornea and/or corneal neovascularization may be indicative of excessive corneal edema.
2. The presence of corneal staining and/or limbal-conjunctival hyperemia can be indicative of an unclean lens, a reaction to solution preservatives, excessive lens wear, and/or a poorly fitting lens.
3. Papillary conjunctival changes may be indicative of an unclean and/or damaged lens.

If any of the above observations are judged abnormal, various professional judgments are necessary to alleviate the problem and restore the eye to optimal conditions. If the CRITERIA OF A WELL FITTED LENS are not satisfied during any follow-up examination, the patient should be re-fitted with a more appropriate lens.

16

6. Adaptation

Visually demanding situations should be avoided during the initial wearing period. A patient may at first experience some mild blurred vision, dizziness, headaches, and a feeling of slight imbalance. You should explain the adaptational symptoms to the patient. These symptoms may last for a brief minute or for several weeks. The longer these symptoms persist, the poorer the prognosis for successful adaptation.

To help in the adaptation process the patient can be advised to first use the lenses in a comfortable familiar environment such as in the home.

Some patients feel that automobile driving performance may not be optimal during the adaptation process. This is particularly true when driving at night. Before driving a motor vehicle, it may be recommended that the patient be a passenger first to make sure that their vision is satisfactory for operating an automobile. During the first several weeks of wear (when adaptation is occurring), it may be advisable for the patient to only drive during optimal driving conditions. After adaptation and success with these activities, the patient should be able to drive under other conditions with caution.

7. Other Suggestions

The success of the monovision technique may be further improved by having your patient follow the suggestions below.

- Having a third contact lens (distance power) to use when critical distance viewing is needed.
- Having a third contact lens (near power) to use when critical near viewing is needed.
- Having supplemental spectacles to wear over the monovision contact lenses for specific visual tasks may improve the success of monovision correction. This is particularly applicable for those patients who cannot meet state licensing requirements with a monovision correction.
- Make use of proper illumination when carrying out visual tasks.

Success in fitting monovision can be improved by the following suggestions:

- Reverse the distance and near eyes if a patient is having trouble adapting.
- Refine the lens powers if there is trouble with adaptation. Accurate lens power is critical for presbyopic patients.
- Emphasize the benefits of the clear near vision in straight ahead and upward gaze with monovision.

* **The decision to fit a patient with a monovision correction is most appropriately left to the eye care professional in conjunction with the patient after carefully considering the patient's needs.**

* **All patients should be supplied with a copy of the Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lens Patient Information Booklet.**

20

PROFESSIONAL FITTING SETS

Lenses must be discarded after single use and must not be used from patient to patient.

WEARING SCHEDULE

The wearing and replacement schedules should be determined by the eye care professional. Regular checkups, as determined by the eye care professional, are extremely important.

Daily Wear

There may be a tendency for the daily wear patient to over-wear the lenses initially. Therefore, the importance of adhering to a proper, initial daily wearing schedule should be stressed to these patients. The wearing schedule should be determined by the eye care professional. The lens is to be prescribed for single-use disposable wear, and is to be discarded after each removal.

MONOVISION FITTING GUIDELINES

1. Patient Selection

a. Monovision Needs Assessment

For a good prognosis the patient should have adequately corrected distance and near visual acuity in each eye. The amblyopic patient or the patient with significant astigmatism (greater than two [2] diopters) in one eye may not be a good candidate for monovision with the Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lenses Contact Lens.

Occupational and environmental visual demands should be considered. If the patient requires critical vision (visual acuity and stereopsis) it should be determined by trial whether this patient can function adequately with monovision. Monovision contact lens wear may not be optimal for such activities as:

1. visually demanding situations such as operating potentially dangerous machinery or performing other potentially hazardous activities; and
2. driving automobiles (e.g., driving at night). Patients who cannot pass their state drivers license requirements with monovision correction should be advised to not drive with this correction, OR may require that additional over-correction be prescribed.

17

HANDLING OF LENSES

When lenses are dispensed, the patient should be provided with appropriate and adequate instructions and warnings for lens care handling. The eye care professional should recommend appropriate and adequate procedures for each individual patient in accordance with the particular lens wearing schedule.

Care for a Sticking (Nonmoving) Lens

If the lens sticks (stops moving), the patient should be instructed to use a lubricating or rewetting solution in their eye. The patient should be instructed to **not** use plain water, or anything other than the recommended solutions. The patient should be instructed to contact the eye care professional if the lens does not begin to move upon blinking after several applications of the solution, and to not attempt to remove the lens except on the advice of the eye care professional.

EMERGENCIES

If chemicals of any kind (household products, gardening solutions, laboratory chemicals, etc.) are splashed into your eyes, you should: FLUSH EYES IMMEDIATELY WITH TAP WATER AND THEN REMOVE LENSES PROMPTLY. CONTACT YOUR EYE CARE PROFESSIONAL OR VISIT A HOSPITAL EMERGENCY ROOM WITHOUT DELAY.

b. Patient Education

All patients do not function equally well with monovision correction. Patients may not perform as well for certain tasks with this correction as they have with bifocal reading glasses. Each patient should understand that monovision, as well as other presbyopic contact lenses, or other alternative, can create a vision compromise that may reduce visual acuity and depth perception for distance and near tasks. During the fitting process it is necessary for the patient to realize the disadvantages as well as the advantages of clear near vision in straight ahead and upward gaze that monovision contact lenses provide.

2. Eye Selection

Generally, the non-dominant eye is corrected for near vision. The following tests for eye dominance can be used.

a. Ocular Preference Determination Methods

- Method 1—Determine which eye is the "sighting dominant eye." Have the patient point to an object at the far end of the room. Cover one eye. If the patient is still pointing directly at the object, the eye being used is the dominant (sighting) eye.
- Method 2—Determine which eye will accept the added power with the least reduction in vision. Place a trial spectacle near add lens in front of one eye and then the other while the distance refractive error correction is in place for both eyes. Determine whether the patient functions best with the near add lens over the right or left eye.

b. Refractive Error Method

For anisometropic corrections, it is generally best to fit the more hyperopic (less myopic) eye for distance and the more myopic (less hyperopic) eye for near.

c. Visual Demands Method

Consider the patient's occupation during the eye selection process to determine the critical vision requirements. If a patient's gaze for near tasks is usually in one direction correct the eye on that side for near.

Example:

A secretary who places copy to the left side of the desk will usually function best with the near lens on the left eye.

18

REPORTING OF ADVERSE REACTIONS

All serious adverse experiences and adverse reactions observed in patients wearing Bausch + Lomb SofLens daily disposable Toric (hilafilcon B) Visibility Tinted Contact Lenses or experienced with the lenses should be reported to:

Bausch & Lomb Incorporated
Rochester, New York 14609
Toll Free Telephone Number
1-800-828-9030
In the Continental U.S., Alaska, Hawaii
1-800-828-9030
In New York State
1-800-462-1720
In Canada
1-888-459-5000

HOW SUPPLIED

Each sterile lens is supplied in a plastic package containing borate buffered saline with poloxamine. Each container is marked with the manufacturing lot number of the lens, diopter power, and expiration date.

21

22